

1 Acid-base titration

EXPERIMENT

Acid-Base Titration

Goal Weak acids are acids that do not dissociate completely, releasing only a few protons into the solution. Acetic acid (CH_3COOH) is a very important weak acid, produced from the fermentation ethanol from the wine. Commercial acetic acid—known as vinegar—is just an aqueous solution of acetic acid. The goal of this experiment is to calculate the *molar concentration* of a sample of acetic acid by means of a standard chemical procedure known as *titration*. In order to do that you will react the weak acid with a basic solution of sodium hydroxide (NaOH), which has a known concentration. You will also use phenolphthalein as an *indicator*.

Background A titration is a technique where a solution of known concentration—often times a base—is used to determine the unknown concentration of another solution—often times an acid. Both substances react with each other in an acid-base reaction. The solution of known concentration is delivered carefully from a buret until an indicator—a third substance added to indicate the end of the titration—changes color. We can express the composition of a solution as the *mass percentage* (%mass) of each component: solute and solvent.

$$\%mass = \frac{\text{mass of solute}}{\text{mass of solution}} \times 100\%$$

For example, if we dissolve 15 g of NaCl in 60. g of water, the total mass of the solution—solute plus solvent—is 75 g and the mass percentage of NaCl in the solution is $(15 \text{ g} / 75 \text{ g}) \cdot 100\% = 20\%$ NaCl . In chemistry, the molar concentration, c , of a solute in a solution, the *molarity* of the solute, is the moles of solute n present in a given volume, V , of the solution in liters. The units of molarity are moles per liter ($\text{mol} \cdot \text{L}^{-1}$), and it is denoted as M .

$$M = \frac{\text{moles of solute}}{\text{L of solution}} = \frac{n}{V}$$

Density, although not a measurement of concentration, is a property of liquids that we can use to relate mass with volume or volume into mass. The formula for density d is:

$$d = \frac{\text{grams of solution}}{\text{mL of solution}}$$

Example

What is the molarity of a sodium hydroxide solution (NaOH , $\text{MW} = 39.997 \text{ g} \cdot \text{mol}^{-1}$) prepared by dissolving 15.00 g of the solute in enough water to make 350.0 mL of solution?

Answer: the molecular mass of NaOH is $39.997 \text{ g} \cdot \text{mol}^{-1}$ and the number of NaOH moles are:

$$n_{\text{solute}} = 15.00 \text{ g} \times \frac{1 \text{ mol}}{39.997 \text{ g}} = 0.3750 \text{ mol}$$

Do not forget to convert the volume from mL to L.

$$350.0 \text{ mL} \times \frac{1 \text{ L}}{1000 \text{ mL}} = 0.3500 \text{ L}$$

The molarity will be:

$$M = \frac{0.3750 \text{ mol}}{0.3500 \text{ L}} = 1.071 \text{ M}$$


Volumetric analysis The determination of concentration by measuring volumes is called *volumetric analysis*. Titrations are volumetric analyses where a buret is used to add and measure the volume of one of the reactants. Acid-base titrations are extensively used chemical techniques employed to determine solute concentration in a solution. In a *acid-base titration*, an acid reacts with a base by gradually adding one solution to the other. The volume of the second solution is known, and the volume of the first solution required for the complete reaction is measured. The formula to use in a titration is:

$$c_a \cdot V_a = c_b \cdot V_b,$$

c_a and V_a are the concentration of the acid and the volume of acid employed, and c_b and V_b are the concentration of the base and the volume of base employed.

A 25 mL solution of perchloric acid, HClO_4 –which has two acidic protons– is titrated with NaOH 0.10 M. The end point for the reaction is reached after 40. mL of the NaOH solution are added. Find the molarity of the acid solution.

Answer: the balanced equation for the acid-base reaction is:

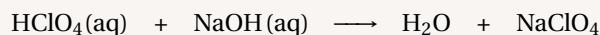


The molarity of the acid, c_a , is unknown whereas the volume used, $V_a = 25 \text{ mL}$, is given.

The base concentration, $c_b = 0.10 \text{ M}$, and volume, $V_b = 40. \text{ mL}$, are given:

$$c_a = \frac{c_b \cdot V_b}{V_a}$$
$$c_a = \frac{0.10 \text{ M} \cdot 40. \text{ mL}}{25 \text{ mL}} = 0.16 \text{ M}$$

An *indicator* is used to indicate the exact end of an acid-base reaction. The indicator chosen for a reaction will have one color before the reaction is complete and a different color when the acid-base reaction finishes. For example in the reaction between acetic acid (CH_3COOH) and sodium hydroxide (NaOH):

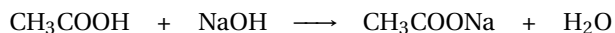

using phenolphthalein as the indicator, the solution will be colorless before completion but pink after completion. At a specific point during the titration, a single drop of the NaOH from the buret will cause the solution being titrated to turn from colorless to a barely discernible pink color. This point is called the *endpoint*.


STUDENT INFO

Name: _____ Date: _____

Pre-lab Done: ☐**Pre-lab Questions**

Acid-Base Titration

1. A 10.00 mL sample of aqueous HCl requires 31.00 mL of 0.0900 M NaOH to reach the endpoint. What is the molar concentration of HCl. The equation for the reaction is:

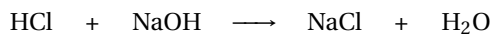

2. The molarity of a vinegar solution is 0.90 M. Calculate the number of acetic acid moles in 10. mL of this solution. Write down your result using scientific notation.
3. Nitric acid (HNO₃) is an acid with three protons. Suppose you titrate 5.00 mL of of this acid with NaOH 0.10 M. Knowing that the end point is reached after 25.00 mL of the base is added, find the molarity of the acid solution.

STUDENT INFO

Name: _____ Date: _____

Pre-lab Done: ☐**Experiment****Acid-Base Titration**

1. Acetic acid titration

- ☐ *Step 1:* – Obtain a 5 mL glass-pipet and a 50 mL buret with a stand and buret clamp.
- ☐ *Step 2:* – Obtain about 30 mL of acetic acid solution (vinegar) in a 50 mL beaker and about 80 mL of the NaOH solution in a clean, dry Erlenmeyer flask. Keep the NaOH solution containing Erlenmeyer closed with a rubber stopper.
- ☐ *Step 3:* – Clean your buret and fill it with the NaOH solution using a plastic funnel.
- ☐ *Step 4:* – Record the initial volume in the buret as 0mL. Read accordingly to the tool precision, including your significant or estimated value.
- ☐ *Step 5:* – Pipet 5.00 mL of acetic acid into a clean 125 mL Erlenmeyer flask that has 20 mL of distilled water and 2 drops of phenolphthalein.
- ☐ *Step 6:* – Record the molarity of the NaOH solution (c_b) indicated in the label of the stock solution bottle. This value will be the same for all experiments.
- ☐ *Step 7:* – Place the flask under the buret. Use a piece of white paper under the flask to distinguish better the color change.
- ☐ *Step 8:* – Add the NaOH solution from the buret in 1 mL portions, while swirling the solution in the flask.
- ☐ *Step 9:* – The titration is completed when an addition of 1 mL causes the color to change from colorless to any shade of pink. Record the final buret volume.
- ☐ *Step 10:* – Repeat the steps above four times and average the resulting acetic acid concentration.
- ☐ *Step 11:* – Make sure you dispose of the solutions and leftovers in the corresponding disposals.

		1	2	3	4
①	Initial Buret Volume (mL)				
②	Final Buret Volume (mL)				
③	NaOH Volume used (mL)				
④	CH ₃ COOH Concentration (M)				
⑤	Mean CH ₃ COOH Concentration (M)				

Calculations

- ① Record the initial volume of the buret. This value is not necessarily 0.00 mL.
- ② Record the final volume of the buret, after you reached the end point.
- ③ The volumen of NaOH used should be: ② — ①
- ④ You can calculate the molarity of the acetic acid solution by means of:

$$c_a = \frac{\textcircled{3} \cdot c_b}{5 \text{ mL}}$$

where c_b is the given molarity of the NaOH solution found in the bottle.

- ⑤ Is the average of the 4 concentrations calculated.

$$\frac{\Sigma \textcircled{4}}{4}$$

STUDENT INFO

Name: _____ Date: _____

Pre-lab Done: ☐**Post-lab Questions****Acid-Base Titration**

1. You need to prepare a sample containing 0.20 g of CuSO_4 from a solution that is 10.% CuSO_4 by mass. What mass of solution do you need?
2. A 10.00 mL sample of aqueous HNO_3 requires 20.00 mL of 0.201 M NaOH to reach the endpoint. Calculate the molarity of HNO_3 .
3. You titrate a vinegar sample—an acetic acid solution in water—with 0.30 M NaOH. Using 10. mL of vinegar, you reach the endpoint after adding 10. mL of the base. Indicate the molarity of the acetic acid solution.